
GERMAN  
POOL

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# neck 頸部按摩器 massager | mn50

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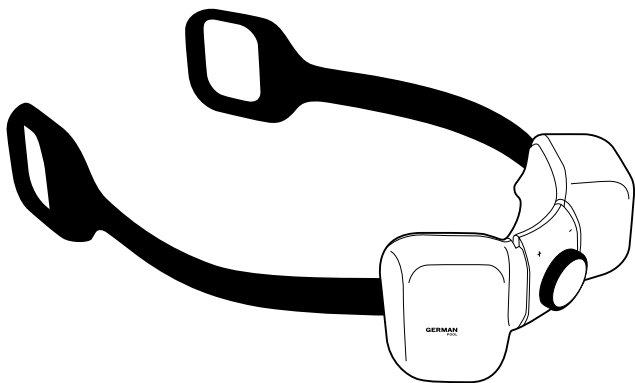---

## IMPORTANT

For safety purpose, you should read the instructions carefully before using the appliance for the first time. Save the instructions for future reference.


The appliance is designed for domestic use and may only be operated in accordance with these instructions for use.

1. Do not use the neck massager while driving.
2. Do not use the massager at a place with high temperature or humidity.
3. Please tie your long hair up before using the neck massager.
4. Please wrap the massager around your neck before switching on.
5. Please use this provided power charger, and pay attention to the safe use of electricity.
6. Do not use this product if you have swellings, inflammation, bruises or any other unusual conditions in your neck. Do not use if you have undergone any neck surgery.
7. Do not use this product if you suffer from hypertension or mental illness.
8. The massager is personal care product, it cannot replace the advice and treatment of health care professionals; in case of any questions or need, please contact your health care professionals directly.
9. This massager is designed for household use, it is not suitable for diagnosis, treatment and other commercial uses.
10. Please wipe the massager with sterile cotton moistened with neutral detergent. Do not clean it with water.


# INSIDE THE PACKAGE

---

This massager comes with a main unit, a power charger and operating instructions.

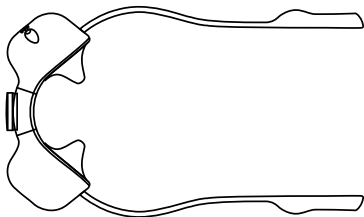Main Unit

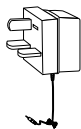Power Charger

Operating Instructions


Protective Sticker

# PRODUCT STRUCTURE

---

## Main Unit
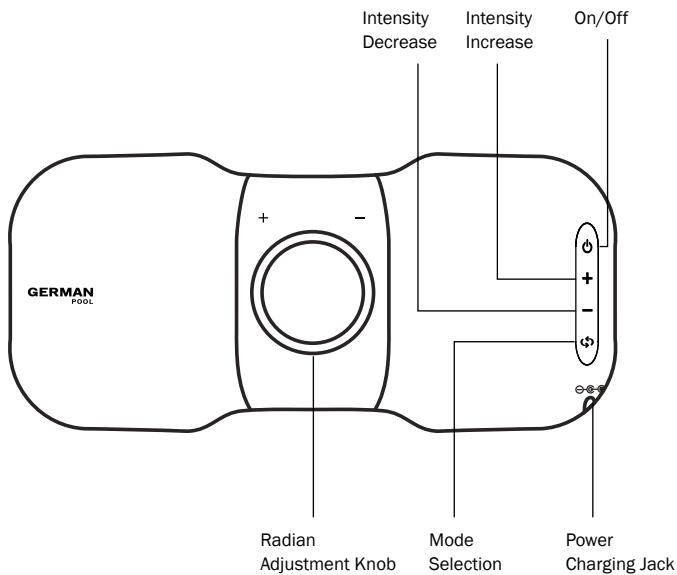
# OPERATING INSTRUCTIONS

---

## Charging

1. Open the plug of the power charging jack at the bottom of the main unit.
2. Plug the DC plug of the power charger into the DC power charging jack of the product.
3. When the product is charging, the LED indicator flashes red/green alternatively; when charging is completed, the LED indicator turns into green.

Note: The main unit is powered by high-performance rechargeable green lithium battery.

## Putting on the Massager

1. Hold the handles with both hands, put the massage head of the main unit on the cervical vertebra, and turn the radian adjustment knob to adjust to the comfortable position for your neck size.
2. Pull the handles with both hands, and adjust the intensity of the massage head contact.

## Switching on/off

### Switching on

Press the  button - LED indicator turns into red, and the product starts operating. (Default: comfortable intensity + shiatsu massage + 20-minute timer)

To switch the massage intensity or massage program, please operate the function keys manually.

### Switching off


Press the  button - LED indicator turns off, and the product stops working (When the massage time is up, the product will automatically switch off)


## Manual mode

Select massage mode

Press the  button - the massage mode changes in sequence: fist percussion massage > kneading massage > pushing massage > shiatsu massage

Select intensity

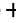Press the  button to switch the massage intensity to slow kneading > comfortable > strong

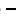Press the  button to switch the massage intensity to strong > comfortable > slow kneading

Note: A beep sound would accompany every touch of a button.

# CARE & MAINTENANCE

---

## Power Check

- Check on the power supply: do not connect the power charger if voltage of the power supply is beyond the prescribed range.

## Maintenance

- Take measures to prevent moisture, rats and pollution.
- Switch off the device before cleaning.
- Put the device in safe place after using. Do not put heavy objects on the device.
- Contact qualified technician if there is any fault. Do not disassemble the device yourself.

# TROUBLESHOOTING

---

If there is any problem when you use the device, please follow the instructions below to check. If the problem still exists or the device cannot work normally, please turn off the power supply and contact service center. Do not disassemble or repair the device yourself.

Do not disassemble or repair the device yourself.

| Phenomenon                    | Possible Cause                                       | Solution                                     |
|-------------------------------|--|--|
| Power-on fail.                | Whether the lithium-battery of the main unit is low. | If yes, please charge the main unit at once. |
|                               | Whether the power charger is not connected properly. | Please ensure it is connected properly.      |
| The massage intensity is low. | If the lithium-battery of the main unit is low.      | If yes, please charge the main unit at once. |

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產品規格及設計如有變更，恕不另行通知。  
Specifications are subject to change without prior notice.

中英文版本如有出入，一概以中文版為準。  
If there is any inconsistency or ambiguity between the English version and the Chinese version, the Chinese version shall prevail.

請填寫以下空格並與發票一併保存，以便維修時核對資料。用戶提供本公司之所有資料必須正確無誤，否則免費保用將會無效。  
Please fill out the form below. This information and the purchase invoice will be required for any repairs. Warranty will be invalid if information provided is found to be inaccurate.

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## Warranty Terms & Conditions

Under normal operations, the product will receive a 1-year full warranty provided by German Pool (effective starting from the date of purchase).

\* This warranty is not valid until Customer's registration information is returned to and received by our Service Centre within 10 days of purchase via one of the following methods:

- 1) Visit our website and register online:  
[www.germanpool.com/warranty](http://www.germanpool.com/warranty).
- 2) Complete all fields on the attached Warranty Card and mail it back to our Customer Service Centre along with a photocopy of the purchase invoice (Warranty Card and invoice copy will not be returned to Customer).

1. Customer who fails to present original purchase invoice will be disqualified for free warranty.
2. Customer should always follow the operating instructions. This warranty does not apply to:
  - labour costs for on-site installation, check-up, repair, replacement of parts and other transportation costs;
  - damages caused by accidents of any kind (including transport and others);
  - operating failures resulted from incorrect voltage, improper operations and unauthorized installations or repairs.
3. This warranty is invalid if:
  - purchase invoice is modified by unauthorized party;
  - product is used for any commercial or industrial applications;
  - product is repaired or modified by unauthorized personnel, or unauthorized parts are installed;
  - serial number is modified, damaged or removed from the product.
4. German Pool will, in its discretion, repair or replace any defective part.
5. This warranty will be void if there is any transfer of ownership from the original purchaser.

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© All rights reserved. Copying, reproducing or using the contents of this manual is not allowed without prior authorization from German Pool, violators will be prosecuted.

[www.germanpool.com](http://www.germanpool.com)